

## RUNWAY® ADJUSTABLE PROSTHETIC FOOT (RS2 or RS2-K1)

*Advanced technology delivers unmatched versatility*

- The only user-adjustable heel height foot designed for active K3 amputees
- Allows the amputee to elevate the heel up to 2" with the push of a button, or by using a 4mm hex key
- Six stiffness categories provide user customization
- Alignment remains consistent, providing exceptionally smooth gait at all heel heights
- Highest weight limit among adjustable heel products—255 lbs
- Ideal product for male and female users seeking versatility and the freedom to wear multiple pairs of shoes of varying heel heights



Runway HX adjustment location



# RUNWAY® (RS2 or RS2-K1)

# RUNWAY® ADJUSTABLE PROSTHETIC FOOT (RS2 or RS2-K1)

**Part #:** Runway (push button) #RS2

Runway HX (hex key) #RS2-K1

**Size:** Standard Toe Shell, **22-28 cm**  
Sandal Toe Shell, **22-28 cm**


**User Weight Rating:** 255 lbs

**Stiffness Category:** 1-6

**Heel Height:** 0-2"

**Connector:** Integral male pyramid

**Warranty:** Carbon Fiber - 36 months  
Ankle Mechanism - 18 months

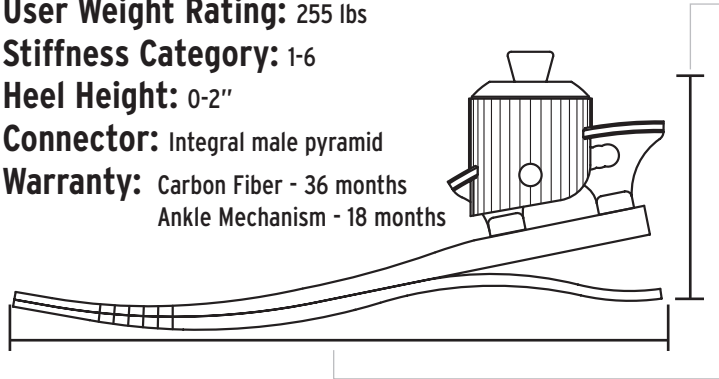

## Clearance & Weight:

Size	Clearance	Runway	Runway HX	Sandal Toe Shell	Mid-Profile Shell
22-25 cm	3.0 in	367g	372g	131g	158g
26-28 cm	3.25 in	449g	414g	210g	210g

Note: A Spectra sock and pyramid connector are included with purchase.

## REQUIRED ACCESSORIES

### Foot Shell Style:

- Sandal-Toe - FTC-2M-S
- Mid profile, standard toe - FTC-2M
- Choice of caucasian or ebony skin tone

## OPTIONAL ACCESSORIES

- Replacement Adjustable Ankle Slide Mechanism- Part #ACC-00-105 (L/R) (available for purchase within the 36 month foot module warranty)
- Foot Shell Removal Tool- Part #AC-00-103
- Additional Shim Adjustment Wrench- Part #AC-00-104
- Additional 4mm hex key for Runway HX adjustments- Part # ACC-00-10400-4M
- Additional Spectra Sock- Part #ACC-00-101xx (8", 10", 12")
- Additional Stiffening Bumpers-
  - KIT-00-119 40 Durometer, Sizes 22-24cm
  - KIT-00-120 40 Durometer, Sizes 25-27cm