

SENATOR® PROSTHETIC FOOT (VS1)

Sophisticated simplicity at an affordable price

- Lightweight, energy-return design facilitates smooth, easy stride and durability
- Foot is not pre-bonded to foot shell, simplifying the use of heel wedges to improve alignment
- Simulated ankle motion provides a stable base of support throughout the stance phase
- Six stiffness categories enhance customization to ensure ultimate comfort and stability
- Ideal for low to moderately active amputees weighing up to 300 lbs (136 kg)



SENATOR® (VS1)

SENATOR® PROSTHETIC FOOT (VS1)

Size*: 22-30 cm

User Weight Rating: 300 lbs


Stiffness Category: 1-6

Connector: Integral male pyramid

Foot Shell Height: Mid-profile
(Included)

Heel Height: 3/8"

Warranty: Carbon - 24 months
Foot Shell - 6 months

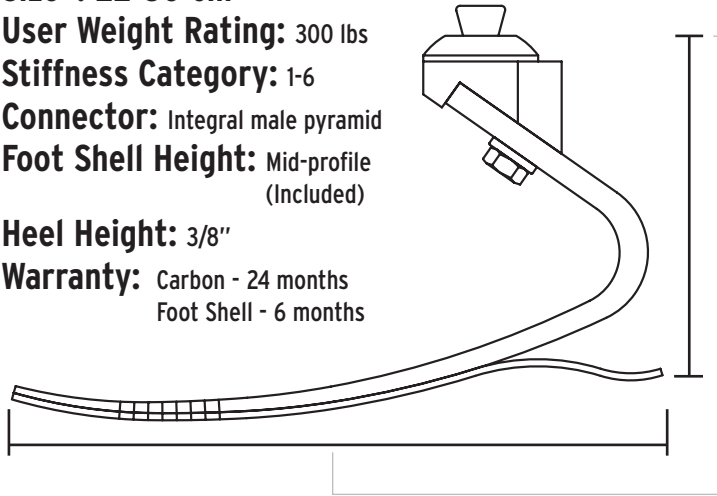

Note: Foot shell, Spectra Sock and connector are included with purchase.

Clearance & Weight:

Size	Clearance	Avg. Weights	
		Foot	Shell
22-25 cm	4.5 in	362g	127g
26-28 cm	5.0 in	390g	159g
29-30 cm	5.4 in	418g	219g

Category Selection Chart:

Weight (lbs)	Impact Level		Weight (kg)
	LOW	MED	
100-120	1	2	45-55
121-145	2	3	56-66
146-175	3	4	67-80
176-210	4	5	81-96
211-250	5	6	97-114
251-300	6		115-136

REPLACEMENT ORDERS & ACCESSORIES

Cosmetic Foot Shell (Replacement): FTC-2M

Choice of caucasian or ebony skin tone

OPTIONAL ACCESSORIES

- **Shell Removal Tool- Part #ACC-00-103**
- **Additional Spectra Sock- Part #ACC-00-101xx (12", 14", 18")**
- **Additional Stiffening Bumpers-**
 - ACC-00-12201 40 Durometer, Sizes 22-25cm
 - ACC-00-12202 40 Durometer, Sizes 26-28cm
 - ACC-00-12203 40 Durometer, Sizes 29-30cm