


SILHOUETTE® PROSTHETIC FOOT (RS4)

High science meets high style

- Slim profile and advanced materials technology facilitate easy and refined cosmesis
- Strong, lightweight carbon-fiber construction allows user to walk farther and faster, using less energy
- Multi-axial rotation delivers an extended range of motion and exceptional stability
- Advanced dual heel system provides a dynamic response for normal ambulation and an override spring for security during higher activity
- Ideal for low to moderate impact K3 amputees
- Rated up to 365 lbs

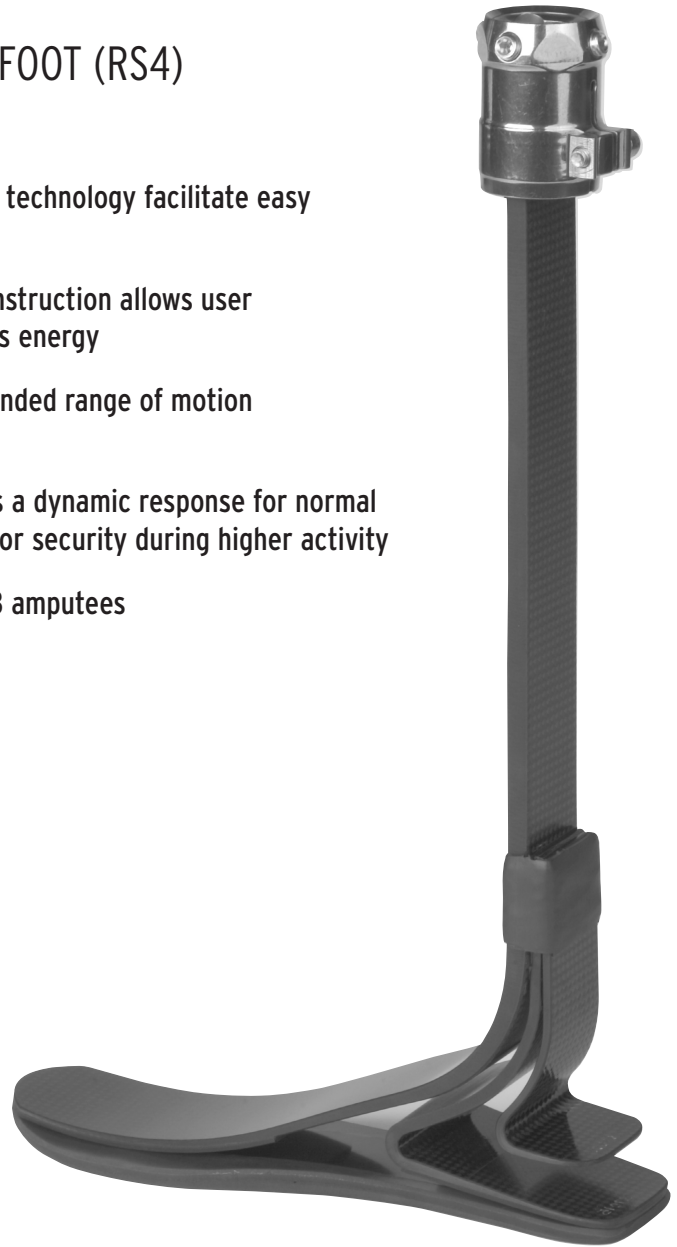

SILHOUETTE® (RS4)

SILHOUETTE® PROSTHETIC FOOT (RS4)

(Incl. 34mm Female Tube Clamp and Insert)

Size*: 22-30 cm

*Special Order Size 31 cm (7-10 day lead time)


User Weight Rating: 365 lbs

Connector: 34 mm female tube clamp with insert

Stiffness Category: 1-9

Heel Height: 3/8"

Warranty: 36 months

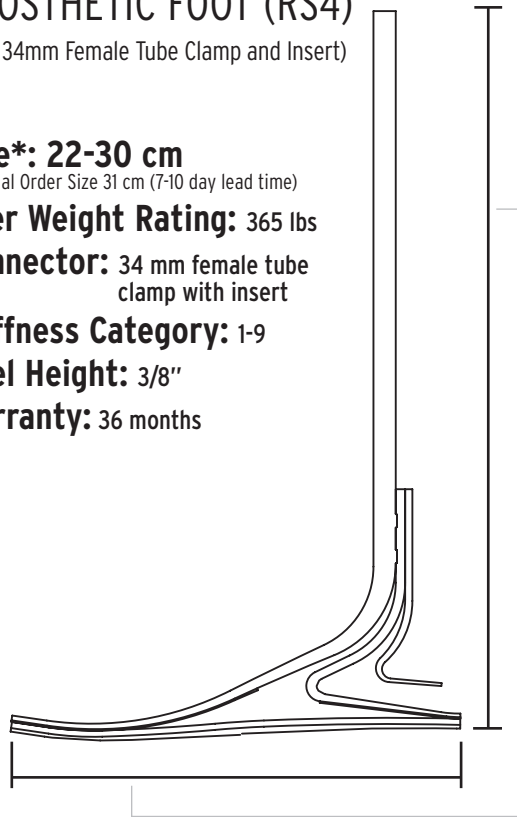

Clearance & Weight:

Size	Minimum Clearance	Avg. Weights	
		Foot	Shell
22-25 cm	7.2 in	506g	127g
26-28 cm	7.7 in	536g	159g
29-31 cm	8.1 in	566g	219g

Note: A Spectra sock and connector are included with purchase.

REQUIRED ACCESSORIES

Cosmetic Foot Shell:

Mid-profile (Default) FTC-2M

High-profile (Optional) FTC-2H

Choice of caucasian or ebony skin tone

OPTIONAL ACCESSORIES

- Replacement 34mm Tube Clamp- Part #ACC-00-123
- Replacement 34mm Female Tube Clamp Insert- Part #ACC-00-124
- Versa Clamp Adapter Kit- Part #KIT-00-123 (includes Extended Tube Clamp and Spacer Rings)
 - Extended 34mm Female Tube Clamp- Part #ACC-EX-123
 - Spacer Ring Set- Part #KIT-00-123S
- Knee Pyramid Adapter, Male (for knees receiving a 30mm pylon)- Part #ACC-00-195
- Replacement Sole Plate, Tri Durometer- Part #RS4-SP
- Foot Shell Removal Tool- Part #ACC-00-103
- Additional Spectra Sock- Part #ACC-00-101xx (12", 14", 18")
- Additional Stiffening Bumpers-

KIT-00-12101 40 & 80 Durometer, Double Sided Tape, & IFU for Silhouette Sizes 22-25cm

KIT-00-12202 40 & 80 Durometer, Double Sided Tape, & IFU for Silhouette Sizes 26-28cm

KIT-00-12203 40 & 80 Durometer, Double Sided Tape, & IFU for Silhouette Sizes 29-31cm