


## SLALOM SKI FOOT (SKI)

*Just step in and go!*

- The first prosthetic foot designed to be quickly and easily inserted directly into a standard ski binding
- Eliminates the need for a ski boot
- Provides necessary flex and shock absorption for fluid motion
- Better control and range of motion increases endurance, allowing more time on the slopes

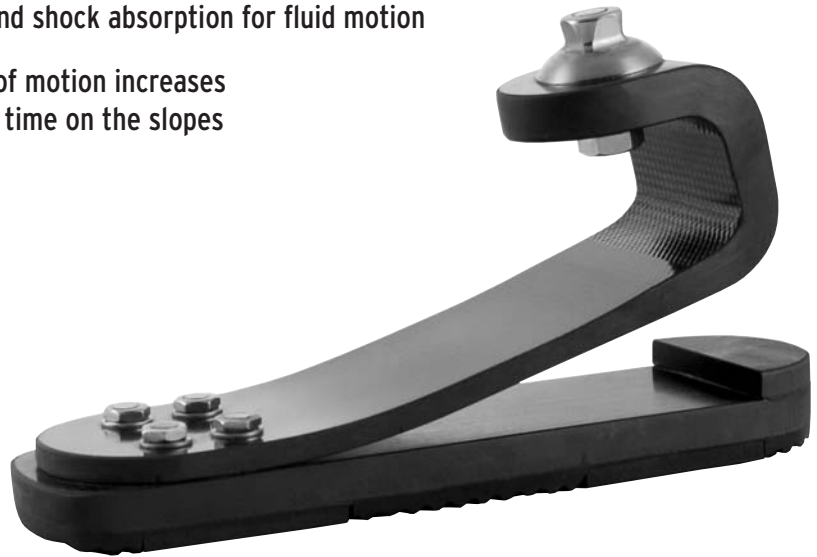

SLALOM SKI FOOT  
(SKI)

# SLALOM SKI FOOT (SKI)

**Size:** Custom made to the size of the ski binding


**User Weight Rating:** 250 lbs

**Stiffness Category:** 11-13

**Connector:** Integral Male Pyramid

**Warranty:** 12 months

**Return Policy:** Non-returnable custom product  
No 60 day trial period

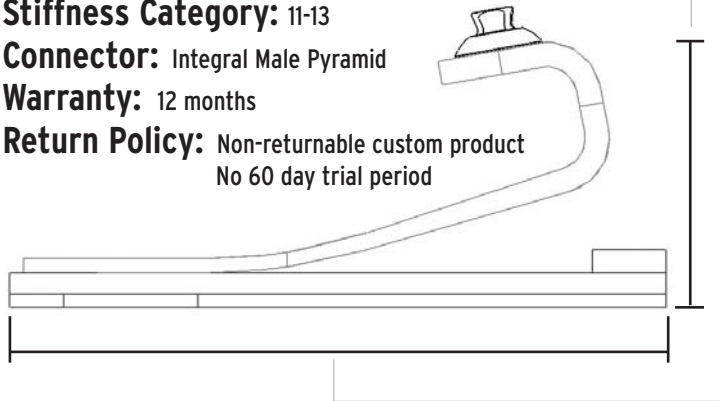

## Clearance & Weight:

Size	Avg. Clearance	Avg. Weight
Custom	5.5 in	1,056g

**Foot Shell Height:** N/A

Please allow 10-14 day lead time.

**Note:** Integral male pyramid is included with purchase.

### Ordering Instructions:

Please provide the following order specifications:

1. Size in cm
2. Patient body weight (Impact level is not required)

### How to measure for size

Unilateral Patient

1. Preferred measurement method
  - a. Measure the length of the patient's sound limb ski boot in cm, from heel to toe.
2. Alternative method
  - a. Measure the length of the ski binding in cm

Bilateral Patient

1. Measure the length of the ski binding in cm

**Note:** The Slalom Ski Foot is made to "DIN specifications" in order to pair with standard ski bindings. If the patient is using something other than standard ski bindings, please contact Customer Service to discuss special needs.

### Category Selection

Weight (lbs)	Category	Weight (kg)
100-160	11	44-72
161-200	12	73-90
201-250	13	91-114